



Face & Body
3D Massage Roller
MINI
User Manual

Distributed By:

Universal Appeal Productions, LLC.
720 South St., #166
Honolulu, HI 96813
www.universalappealproductions.com
www.3Dmassageroller.com



www.3Dmassageroller.com

STRUCTURE DESCRIPTION

- The unique V-shaped structured design with 360° multi-angle rotating kneading technology has a good sense of grip and deep squeeze to lift the skin.
- Diamond-shaped 3D rolling balls can be fit to all parts of body and facial curves.
- Micro Current solar panel on the handle, absorbs energy from the light, transfers a low electric current which flows to your body during massaging.
- Ergonomically streamlined massage handle can be used to easily accommodate working on different massage spots.
- Waterproof, simple to use, portable, no batteries or charging cord required, easy to clean.

BENEFITS

- Promotes face lift, skin tightening and toning, body shaping and slimming.
- Helps to diminish the appearance of puffiness, fine lines, wrinkles.
- Improves blood circulation, alleviate tension.
- Exercises muscle tissue and loosens tight spots, releases stress, relieves fatigue.
- Enhances skin elasticity and muscle firmness to prevent sagging.
- Raised triangle texture can effectively stimulate accupuncture points to help relax.
- Hugs and squeezes every curve of your face and body.

TIPS FOR BASIC USAGE


- Use 5 minutes per treatment area per day.
- Increase or decrease pressure by lowering or raising the end of handle.
- 3D Massage roller can be used with oil, lotion or skin cream for massage and helps to absorb.
- Can be cleaned with water or cloth.
- You can use 3D Massage Roller on chin, cheeks, forehead, neck & shoulder, arms & fingers, back, waist, belly, buttocks and thighs, legs.

FACIAL CARE

- For cheeks: Rolling along the facial line from bottom up to the mouth. Tighten the corners of the mouth, cheeks and facial skin.
- For chin: Focus on massaging chin can remove double chin effectively. Don't overexert.
- For forehead: Rolling horizontally from one side to other side, or vertically above the eyes to head, releases the stress and alleviates tension.
- For neck: Grip with the handle up, scroll up and down, back and forth from the shoulder along the back of the neck.

BODY CARE

- For arms: Roll back and forth from elbow to armpit to reduce the flappy muscles (reduce/remove redundant fat).
- For waistline: Scroll up and down the sides of the waist along the lateral lumbar curve.
- For chest: Scroll to the center of the chest from the bottom of the chest and from armpit, good for breast enhancement and shaping.
- For legs: Roll back and forth in the inner side and the back side of thighs.

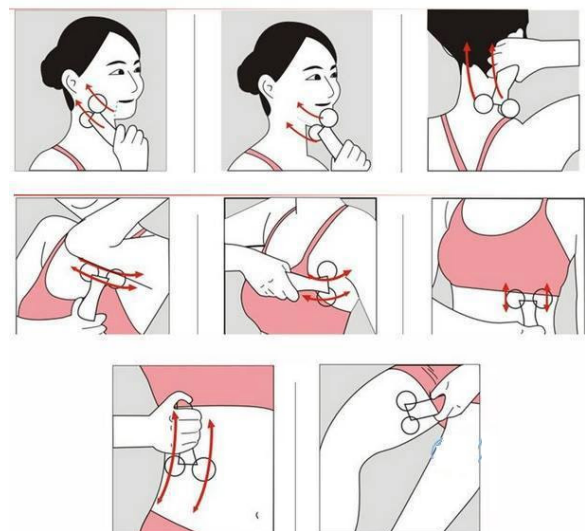

CAUTIONS

- Please remove necklace, earrings, and hand accessories before use to prevent getting tangled with the device.
- Please tie up hair, keep it fixed when in use to avoid hair curling in the device.
- Please do not use on protruding parts like spinal bones.
- Please do not use with too much strength or overexert.
- Please keep out of children's reach.
- Please do not use if you are allergic to metal.

DISCLAIMER

- Use this device at your own risk.
- Neither Universal Appeal Productions, LLC. or any other distributors and retailers assume responsibility or liability for any injuries or damages resulting from using this device.

PRODUCT SPECIFICATIONS

Face & Body 3D Massage Roller Mini
Size: 100 X 54 X 38 cm
Weight: 135 g
Material: ABS + Metal



Face & Body
3D Massage Roller
User Manual

Distributed By:

Universal Appeal Productions, LLC.
720 South St., #166
Honolulu, HI 96813
www.universalappealproductions.com
www.3Dmassageroller.com



www.3Dmassageroller.com

STRUCTURE DESCRIPTION

- The unique V-shaped structured design with 360° multi-angle rotating kneading technology has a good sense of grip and deep squeeze to lift the skin.
- Diamond-shaped 3D rolling balls can be fit to all parts of body and facial curves.
- Micro Current solar panel on the handle, absorbs energy from the light, transfers a low electric current which flows to your body during massaging.
- Ergonomically streamlined massage handle can be used to easily accommodate working on different massage spots.
- Waterproof, simple to use, portable, no batteries or charging cord required, easy to clean.

BENEFITS

- Promotes face lift, skin tightening and toning, body shaping and slimming.
- Helps to diminish the appearance of puffiness, fine lines, wrinkles.
- Improves blood circulation, alleviate tension.
- Exercises muscle tissue and loosens tight spots, releases stress, relieves fatigue.
- Enhances skin elasticity and muscle firmness to prevent sagging.
- Raised triangle texture can effectively stimulate accupuncture points to help relax.
- Hugs and squeezes every curve of your face and body.

TIPS FOR BASIC USAGE


- Use 5 minutes per treatment area per day.
- Increase or decrease pressure by lowering or raising the end of handle.
- 3D Massage roller can be used with oil, lotion or skin cream for massage and helps to absorb.
- Can be cleaned with water or cloth.
- You can use 3D Massage Roller on chin, cheeks, forehead, neck & shoulder, arms & fingers, back, waist, belly, buttocks and thighs, legs.

FACIAL CARE

- For cheeks: Rolling along the facial line from bottom up to the mouth. Tighten the corners of the mouth, cheeks and facial skin.
- For chin: Focus on massaging chin can remove double chin effectively. Don't overexert.
- For forehead: Rolling horizontally from one side to other side, or vertically above the eyes to head, releases the stress and alleviates tension.
- For neck: Grip with the handle up, scroll up and down, back and forth from the shoulder along the back of the neck.

BODY CARE

- For arms: Roll back and forth from elbow to armpit to reduce the flappy muscles (reduce/remove redundant fat).
- For waistline: Scroll up and down the sides of the waist along the lateral lumbar curve.
- For chest: Scroll to the center of the chest from the bottom of the chest and from armpit, good for breast enhancement and shaping.
- For legs: Roll back and forth in the inner side and the back side of thighs.

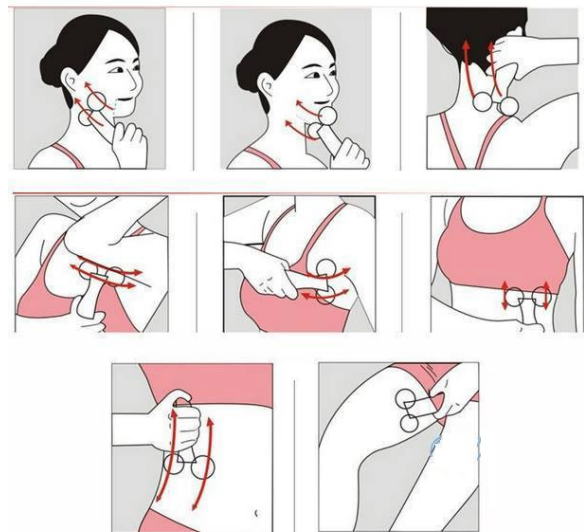

CAUTIONS

- Please remove necklace, earrings, and hand accessories before use to prevent getting tangled with the device.
- Please tie up hair, keep it fixed when in use to avoid hair curling in the device.
- Please do not use on protruding parts like spinal bones.
- Please do not use with too much strength or overexert.
- Please keep out of children's reach.
- Please do not use if you are allergic to metal.

DISCLAIMER

- Use this device at your own risk.
- Neither Universal Appeal Productions, LLC. or any other distributors and retailers assume responsibility or liability for any injuries or damages resulting from using this device.

PRODUCT SPECIFICATIONS

Face & Body 3D Massage Roller
Size: 15.5 X 9 X 5.30 cm
Weight: 168 g
Material: ABS